



TheraPro

SERIES

BEST
PRACTICES

Full Spectrum InfraRed
Applied Sciences



Table of Contents

Thank you and Welcome!	1
Congratulations	2
What's In The Box	3
Setup Guide	4
How to Use	5
How to Get Started	6
Benefits & Potential Reactions	7
Core Technologies	8
Remote	9
PEMF Levels	10

How to Use PEMF	11
How to Use TENS	12
Travel Instructions	13
Special Populations	14
Specs	15
Specs	16
Celsius to Fahrenheit Chart	17
FAQ's	18
Warranty	19
Contact	20

Clean Manufacture

Our materials have been specially selected and tested to assure they are free of any toxins and 3rd-party tested to meet the highest international standards for bio-compatibility.

Very Low EMFs/ELF/RF

You've got enough EMFs to worry about. All of our technology is developed with as low EMF rating as possible.


Thank you & Welcome!



Whether this is your first Therasage product or the next addition to your health collection, thank you for choosing us in your journey to better health and well-being!

Therasage began from our personal journey to help our daughter manage, and eventually overcome her symptoms from Chronic Lyme Disease. With all of the challenges of our current healthcare approaches, it is our vision that each of us takes more responsibility for our health needs by incorporating a positive, natural, healthful daily regimen. We can make a positive statement toward better health using healing technology.

Therasage's mission is to deliver the best and the most advanced devices and applications in the natural personal health space. We hope that many of our FDA-listed rating, coupled with our commitment to corporate responsibility, integrity, and the environment makes us a worthy partner in your journey toward optimal health. Welcome!

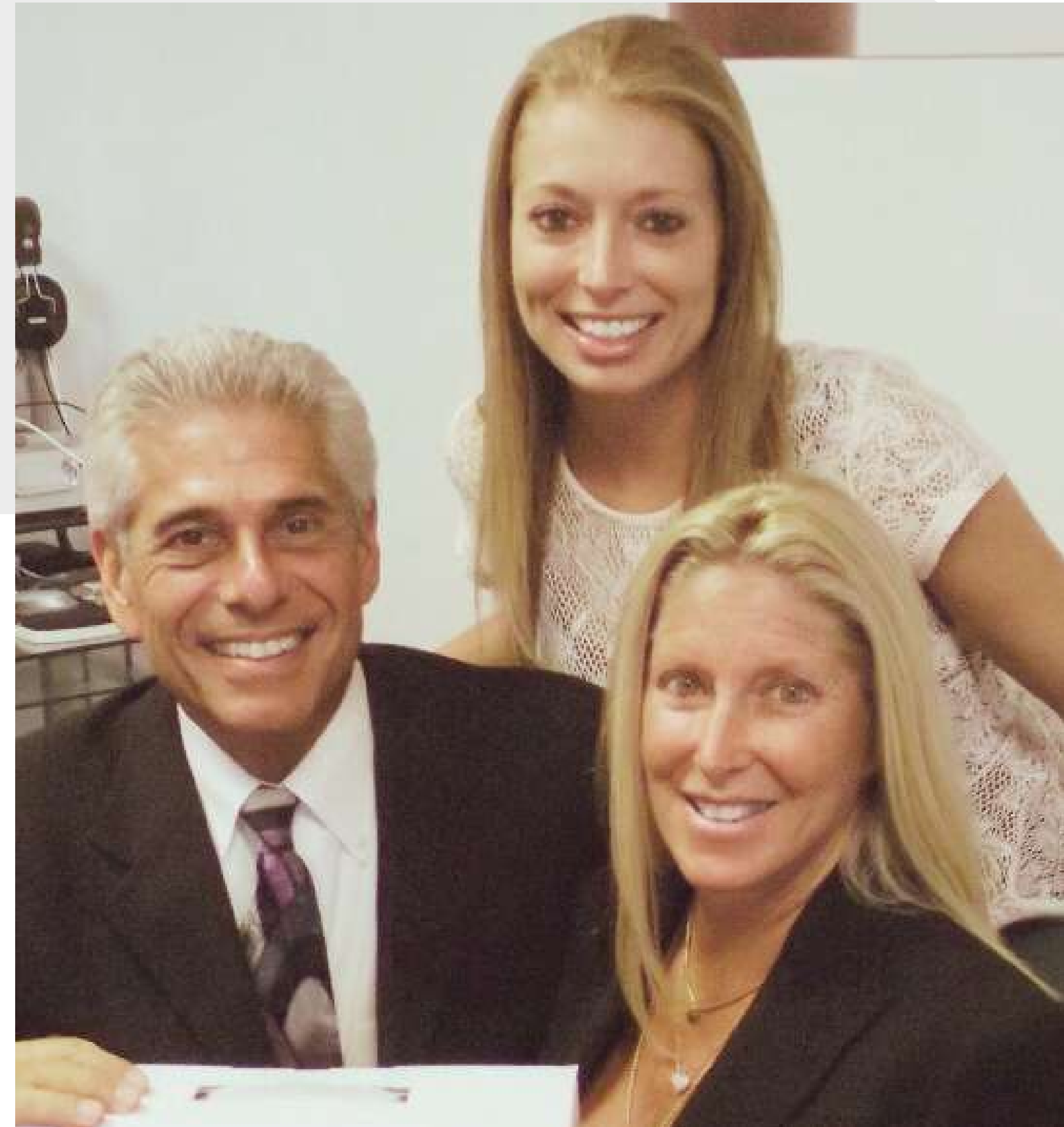

TheraPro

Congratulations on your purchase of the TheraPro Healing Pad!

At Therasage, we pride ourselves as pioneers in advanced wellness technologies. Our commitment to innovation has led us to develop the TheraPro Healing Pad, a unique device that integrates six powerful healing modalities into one comprehensive treatment solution.

The TheraPro Healing Pad is engineered to revolutionize your health regimen by providing PEMF, full-spectrum infrared, gemstone therapy, advanced earthing, red light therapy, TENS for pain relief, EMF remediation, and bio-energetic harmonizing technologies. This unparalleled combination ensures that you can harness the benefits of multiple therapies simultaneously, a feat previously not achievable.

By incorporating the TheraPro Healing Pad into your daily wellness routine, you can enjoy enhanced cellular function, increased circulation, reduced inflammation, and accelerated recovery from physical discomfort. Research supports the synergistic effects of these technologies in promoting overall health and vitality.

Embrace the TheraPro Healing Pad as a cornerstone of your health practice and experience the transformative effects of Therasage's innovation. Here's to a healthier, more vibrant you!

"I use my TheraPro everyday, throughout the day. I have it on my chair while I'm working, or setup next to my couch when I'm relaxing watching TV. I find the more accessible it is the more I use it."
- Joline

"PEMF therapy does not care what YOU perceive to be wrong w/ your body, it can provide support to you whether you have a broken bone, failing heart, anxiety or any other health ailment. It can help stabilize the body to address underlying imbalances to almost any health condition."
- Ben

This reference guide is designed to help you get the most out of your TheraPro Healing Pad. Our approach ensures a gentle transition for your body, creating a healthy foundation for long-term wellness.

****We are NOT offering medical advice**--these are merely guidelines.





TheraPro

Your package includes:

- 1 x TheraPro Healing Mat (in either large or regular)
- 1 x Canvas Carry Bag
- 1 x TheraPro Remote
- 1 x 6ft TENS Cord
- 2 x TENS Pads
- 1 x User Manual

Please check the contents to ensure everything is included.

Step 1: Unpack your TheraPro Healing Mat and make sure all parts are accounted for. If you are missing a part, please email us at support@therasage.com or call us at 1-888-416-4441.

Step 2: Read these pages:

- How to Use: MUST READ
- How To Use TENS
- How To Get Started

Step 3: **Connect the Device:** Plug your cord into the mat as shown.



Step 4: For the first few uses keep the remote pages nearby! Soon you will become familiar with the buttons, and the PEMF levels.

Step 5: Turn the **Power On** by pressing the silver circle, select and set your programs per the **Remote** and the **PEMF Levels** pages.

Step 6: Relax and Enjoy Your Session!

How to Use: MUST READ



Welcome to a transformative experience with your TheraPro Healing Pad! Follow these steps to maximize the benefits and enjoy a seamless therapy session:

Ready When You Are

Keep the Healing Pad folded and close by—whether next to your favorite couch or chair, it's designed for convenience. Always keep it plugged in during use to maintain your settings, but be sure to unplug it when not in use to conserve energy.

Direct Contact

Place the gemstone side against your body. If you're wearing clothes, that's fine! However, avoid wearing protective or thick clothing as it can filter and reduce the effectiveness of the therapy. Materials like latex or other stretchy fabrics may diminish the therapeutic results.

Hydrate! Hydrate! Hydrate!

Hydration is crucial! Make sure to drink plenty of water within one hour before and after your session. The infrared therapy mobilizes toxins, making hydration essential to support the detox process effectively. Do not underestimate the importance of this step—the TheraPro is a potent detox tool.

All Together Now

Utilize all available therapies at once for maximum benefit, except for TENS, which should be used more cautiously and specifically. Check the TENS precautions to ensure it's suitable for your needs before starting.

Use It Daily, But Start Slow

Incorporate the Healing Pad into your daily routine to leverage the cumulative benefits of the therapies. Start with shorter sessions and gradually increase your usage as recommended on the Getting Started page, allowing your body to acclimate comfortably.

It's Rude To Stare...And Not Recommended

Though it's safe, avoid staring directly into the red lights. Like looking at the sun, it's not harmful in short bursts, but it's not polite—or wise—to stare.

How to Get Started



We know you want to jump right in and start enjoying the benefits of this powerful healing tool...we just want to make sure you start slow and acclimate yourself so you can avoid any uncomfortable symptoms of detox (can we say HYDRATE again?)

Keep It Low & Slow For 2 Weeks & Hydrate!

ESSENTIAL - Do NOT skip the acclimation period!

80% will be able to add more sessions after 2 weeks. Do what feels right for you after 2 weeks. If using less than twice per day, it may take longer to acclimate. *Check the **Precautions** and **Special Population** pages before starting and for special instructions.

WEEKS 1-2

Use up to twice per day, for 1 hour at a time (other than TENS - see TENS instructions + precautions pages)

1

WEEKS 3-4

Use all the technologies as often as you'd like (other than TENS - see TENS instruction page)

2

WEEKS 5+

Give your body Thera**Love** as much as you'd like! Play with the different PEMF settings (see PEMF Levels page) for the potential benefits at each level

3

D R I N K W A T E R ! !

Benefits & Potential Reactions

Research has shown that the healing technologies incorporated in our TheraPro Healing Pad offer numerous benefits that support the body's natural healing and regenerative processes. Regular use can provide soothing warmth, reduce pain, enhance relaxation, boost circulation, and alleviate stress and anxiety. It can also amplify the effects of other integrative therapies. The potential benefits are extensive and the positive effects accumulate over time—so start early and use often!

Each technology delivers unique advantages, and when used together, their combined effects can be transformative. The benefits at the cellular level and improvements in sleep quality can positively impact every organ and system in your body.

You Might Enjoy:

- Short and long-term pain relief
- Improved circulation
- Lowered stress and anxiety
- Improved relaxation
- Healthier skin
- Improved detoxification
- Increased sports performance
- Faster recovery times
- Higher levels of oxygen in your blood
- Reduced muscle tension
- Improved sleep
- More energy
- Clearer thinking
- Improved mood
- A more balanced system
- Improved brain function
- Enhanced nutritional absorption
- Increased blood oxygenation



Some Might Experience Symptoms of Detox:

As your body acclimates and circulation increases, your cells, and organs release toxic buildup. This can potentially cause detox reactions that could include:

- Headaches
- Changes in bowel movements
- Fatigue
- Nausea
- Sluggishness
- Skin breakouts
- Moodiness

Important Note: Detox reactions are different for everyone.

What Can You Do?

- Add organic lemon/lime to your water
- Jump in your Infrared Sauna to help sweat out the toxins
- Take a hot bath with Epsom salt

HYDRATE! HYDRATE! HYDRATE!
(We can't emphasize this enough)

Core Technologies



PEMF Therapy (Pulsed Electromagnetic Field)

PEMF has been around for 60+ years and is a rehabilitation therapy that helps regenerate damaged or injured cells, increase microcirculation, and support the energy factories within each cell (the mitochondria) to re-energize & give your body what it needs to heal on a cellular level. Nine program choices within PEMF give a huge range of potential benefits.

Full Spectrum Infrared

The core technology also in our Infrared saunas, InfraRed heat penetrates deep into your body to support improved circulation, pain relief, detoxification, stress relief, improved immune function, and much more!

Red Light

Uses low-level wavelengths of red light to stimulate skin cells and promote healing.

Gem Stone Technology

We use Natural Rose Quartz, Jade Stone, Tourmaline, Shungite and Amethyst that are all Piezoelectric Earth Materials. They have been used for thousands of years to provide negative ions and restored energy fields to the body.

TENS (Transcutaneous Electrical Nerve Stimulation)

TENS can also be considered a Micro-current Electron donor device, sending electron pulses through the skin to charge up your cell battery and kickstart your body's own painkillers. Think of E-Stim devices you've used at Physical Therapy office.

Earthing/Grounding

Earthing equalizes your body's charge with the electric potential of the earth and improves the voltage level in your body. This stability counteracts the fluctuating electronic fields around us, supports a strong immune system and improves cell communication.

EMF Remediation – Harmonizing – Shielding & Blocking

Our bio-energetic technology has two levels of support. The first level of action is direct—it harmonizes the negative energetic frequencies of electromagnetic radiation to minimize its health impact. The second level is adaptive—it supports the cells that are stressed by cell tower radiation, bringing their energy to an optimal level for good health and well-being.

Remote

Familiarize yourself with the remote—it might seem like there are many buttons at first! Keep the next two pages handy until you feel comfortable navigating through the programming.

- 1 TENS Plug In Socket:** For TENS pad connection.
- 2 Heat Intensity Control:** Press (+) or (-) to adjust the temperature. Press and hold together for 1 second to switch between °F and °C.
- 3 Timer Control:** Press to select your time preference between 1 and 12 hours. With each press, it will reduce the time by one hour.
- 4 Red Light Control:** Press to turn on the Red Light Therapy function. The red light will light up on the pad and will run for 2 hrs.
- 5 PEMF On/Off/Selection:** Press once to turn on. Click through programs P1-P9. See next page for specifics of each level – Turns off after 30 minutes.
- 6 Power On/Off:** Press and hold for 3 seconds to power on. The system will automatically detect harmful EMF (Electromagnetic Waves Field) and eliminate them. Press once to turn off.
- 7 TENS On/Off and Strength Adjustment:** Press once to turn on. You can modify the intensity when blinking.
- 8 TENS Intensity Control:** Press the **TENS INTENSITY** (+) or (-) buttons to adjust intensity level.
- 9 Time/Settings Display:** Shows time setting for TENS or PEMF treatment, and PEMF level. See activated symbols above display.
- 10 Temperature Display:** From 50°F – 158°F degrees or 10°C – 70°C. Use the **HEAT INTENSITY** (+) and (-) buttons to adjust. Takes time to heat up. Look at a Celsius-Fahrenheit chart for more accurate conversions.



PEMF Levels

& Potential Benefits Of Each Level

1 PEMF On/Off and Program Setting:

- Press for 1 second to turn on. Works for 30 minutes, then automatically turns off.
- Click this button repeatedly to move through the program levels. Each click will light up the corresponding program identifier to the left.
- The selected program will be confirmed within 3 seconds, and will show in the LED screen on the top right.

2 PEMF Program Details

in order of appearance, top to bottom

R-P9 : Randomized mix of all frequencies; cycles through programs 1-8, with 6-second interval change of frequency (shown as **R** in illustration).

64Hz-P8 :
56Hz-P7 : } Pain Relief
48Hz-P6 :

40Hz-P5 : (Gamma Waves, γ Waves) Peak concentration, brain optimization

25Hz-P4 : (Beta Brain Waves, β Waves) Energize

10Hz-P3 : (Alpha Waves, α Waves) Calm

8Hz-P2 (7.83Hz) : (Schumann Resonance)
 Mimics earth's magnetic energy exactly, Akin to grounding barefoot to the earth.

3Hz-P1 : (Delta Waves) Repair/Restore Mode
 (shown as circled **3** in illustration)

Potential Benefits

P1 Improves deep sleep.

P2 Healing, blood circulation and lowers inflammation.

P3 Lowers Stress, reduces anxiety and depression and increases creative thinking

P4 Alertness

P5 Improves Memory

P6 Level 1 – Healing, blood circulation and pain relief

P7 Level 2 – Healing, blood circulation and pain relief

P8 Level 3 – Healing, blood circulation and pain relief

P9 Mix of all benefits mentioned above



THERASAGE®

How To Use PEMF

PEMF therapy can be felt differently by each person. While some may not feel anything, others might experience a slight warmth in areas with increased blood flow and enhanced microcirculation. Don't be surprised if you notice benefits both during and shortly after your session begins—some people feel a surge of energy, while others may feel sleepy. The response varies from person to person. Let's dive in and learn how to use this amazing technology!

First Of All...CAN You Use It?

PEMF is not recommended if you are pregnant or have given birth within the last 2–4 weeks. For other specific conditions, please refer to our “Special Populations” section on page 15. We advise consulting with your healthcare provider or contacting our Thera Optimal Health team at support@theraoptimalhealth.com for guidance before starting this therapy.

Ok, Now Let's Turn It On!

To activate, press button #1 (shown on previous page) on your device for 3 seconds. Continue pressing to select your desired intensity level. promote healing.

Choosing Your Level

Select the treatment level based on your goals. Use 3Hz, P1 for better sleep or 25Hz, P4 for increased alertness. For a full list of levels and their benefits, refer to the previous page.



Recommended Usage

For beginners, start with one hour sessions twice per day for the first two weeks. After acclimating, adjust the duration and frequency as desired.

Where To Use On The Body

Apply the PEMF device to any area of your body that needs Thera**LOVE**. It's designed to be moved around to target various spots that could benefit from increased blood flow.

Important Precautions

HYDRATE!! Failure to drink enough water could lead to detox symptoms. Also, consult the “Special Populations” section on page 15 to ensure PEMF therapy is suitable for you at this time.

How To Use TENS:

TENS (Transcutaneous Electrical Nerve Stimulation) is a potent modality for pain relief, which can be used alongside other TheraPro modalities. However, it comes with specific precautions and contraindications that you **must read carefully** to ensure it's safe and beneficial for you.

Who Should NOT Use TENS:

Avoid using TENS if you have any implanted metallic or electronic device, such as a spinal cord stimulator, cardiac pacemaker, defibrillator, or any in-dwelling pumps or monitors. Consult your healthcare provider if you have conditions like cancer, diabetes, epilepsy, cognitive impairments, or if you are pregnant.

Precautions: How NOT to Use TENS

- Do not apply TENS to the head, face, throat, swollen, red, infected or inflamed skin, cancerous lesions, open wounds, rashes, or areas lacking normal sensation.
- Avoid using on your spine directly or simultaneously on the front and back of your chest.
- Never use TENS while bathing, showering, or sleeping.
- Keep out of reach of children.
- Discontinue use and consult your healthcare provider if you experience skin irritation, dizziness, nausea, or headaches.

OK, now that you CAN use TENS:

Connect the cord to the TENS socket on your remote. Peel off the sticky backing of the application pads and apply directly to clean, dry skin, avoiding any lotions or oils.

Turn It On!

Turn on the device by pressing the TENS On/Off button (see button #7 on page 9). Adjust intensity using the (+) and (-) buttons. Start at a low intensity and increase gradually to ensure comfort.

Pad Placement

Position the TENS pads on either side of the painful area but never touching each other or any metal object you might be wearing. They should be spaced about **3-6 inches apart**.

How Long

Begin with a **30-minute** session at the lowest intensity and adjust slowly to a comfortable level.

Adjustment

The sensation should be a gentle tingling; ensure it remains comfortable at all settings. As you become accustomed to TENS, you may require higher intensity to achieve the same level of stimulation. Always ensure it remains comfortable.

Maintenance

After use, turn off the TENS by pressing the TENS on/off button until the light and function turn off.

Travel Instructions



Healing to go! An active life often includes travel—whether for work or while on holiday you don't need to interrupt your healing and calming sessions. There aren't many attachments, the controls and connectors are durable, and the gemstones and wiring are protected by folding the pad in on itself.

Cool it down! Wait until the pad is cooled down after use to wrap and pack.

The **Pro-Regular** Pad is the easiest to take with you when flying. It folds into a manageable size and slips into the provided bag. You can easily carry it on or slide it between the clothing in a checked bag. It DOES fit into a rolling carry-on bag.

The Pro-Large is just as portable, and though it will fit into some checked luggage, size and weight are significant. A road trip may be this unit's best means of travel. Strap it onto the car seat, the controller can also be wall-mounted.

➔ International Plugs

Your TheraPro Series Healing Pads use 110V, and will work with any international converter that allows plugging into a 110V wall socket.

The TheraPro is also available in 220V, specifically for use in countries that use 220V.

➔ Size

The **Pro-Regular** is 33.5" x 18.5" | 85cm x 48cm when laid out. Folded dimensions are 18.5" x 12" x 3" | 47cm x 30cm x 7.6cm

The **Pro-Large** is 58.5" x 26" | 148cm x 66cm when laid out. Folded dimensions are 18.5" x 12" x 3" | 47cm x 30cm x 7.6cm

➔ Weight

The **Pro-Regular** weighs 9 lbs or 4.08 kg.
The **Pro Large** weighs 24 lbs or 10.88 kg.



The TheraPro Healing Pad can be used by almost anyone, however please be aware that Special Populations will have specific precautions and limitations depending on the modality.

Make sure to follow acclimation guidelines on the How To Get Started Page.



Fragile Persons of Any Age

If you are very symptomatic, or have chronic challenges such as fatigue, Lyme, long COVID, etc., please always contact your primary healthcare provider before using any new technology. Also, contact our Thera Optimal Health experts to explore further how you can benefit from this product. Contact us at support@theraoptimalhealth.com

Anyone Under 18 Years Of Age

We cannot decide what is best for your child. We can provide basic guidelines, but ultimately it is up to you and your healthcare provider to decide.

Infants and Children 0–6 years

- Can use all modalities (except TENS) for a 15-minute treatment with an adult present.
- TENS: Please check with Thera Optimal Health staff first as to why you would want to use it, at support@theraoptimalhealth.com.

Kids 6–18 years

- Generally, you can treat with half the time as an adult (15–30 minutes). For more specific use contact Thera Optimal Health (support@theraoptimalhealth.com).

Pregnant People

- Do NOT use PEMF or TENS while pregnant and for 2–4 weeks after delivery.
- Red and Infrared Light have proven to be effective--there is no science suggesting that biophoton modulation is contraindicated for pregnancy--but it is still not recommended.

This reference guide is meant as a starting point to enhance your TheraPro experience. Our approach is conservative to allow your body to gently acclimate and create a healthy natural platform for long-term good health and well-being.

We are NOT offering medical advice--these are merely guidelines. There are more manicured approaches based on special needs. For more specific guidance, contact us to make an appointment with our Thera Optimal Health advisory team and always check with your healthcare provider.

Category	TheraPro Regular	TheraPro Large
PEMF	9 – Settings – 3-64 HZ	9 – Settings – 3-64 HZ
IR	Full Spectrum – NIR – Mid – FAR – 700nm-14um	Full Spectrum – NIR – Mid – FAR – 700nm-14um
Red Light	YES – 663nm / 2hrs Auto Shut off	YES – 663nm / 2hrs Auto Shut off
Neg Ions	YES	YES
Gem Stones	Shungite, Jade, Amethyst, Tourmaline, Rose Quartz	Shungite, Jade, Amethyst, Tourmaline, Rose Quartz
Grounding	Yes – Electrical	Yes – Electrical
TENS	YES	YES
Controller	Digital – See Below Details – Remembers Last Settings	Digital – See Below Details – Remembers Last Settings
Wave Form	Square & Sine	Square & Sine
EMF	Low	Low
Temperature	86 F – 158 F / 10 C – 70 C	86 F – 158 F / 10 C – 70 C
Time	1-12 hours	1-12 hours
# of Coils	4	6
Pulse Duration	3 Drop Off Rates	3 Drop Off Rates
Materials	NON Toxic – No VOC's	NON Toxic – No VOC's
Size	18 x 33	26 x 58
Strapping	YES	YES

Category	TheraPro Regular	TheraPro Large
Flexible	YES	YES
Electrical	110v / 220v	110v / 220v
Carry Case	YES	YES
Power Rating	966 gauss	1100 gauss
Weight	9 lbs	24 lbs
TENS	90v – 110v – 120min auto shut off – 20 intensities	90v – 110v – 120min auto shut off – 20 intensities
TENS Pads	Reusable / Accessory	Reusable / Accessory
PEMF	30 mins Auto Shut Off	30 mins Auto Shut Off
	P1 – Delta – 3hz	Parasympathetic
	P2– 7.83hz	Schuman Resonance – grounding – lowers inflammation
	P3 – 10hz	Alpha waves – lowers stress & anxiety
	P4 – 25hz	Beta waves – sympathetic – alertness
	P5 – 40hz	Gamma waves – memory
	P6 – 48hz	Pain Relief level 1 – blood circulation
	P7 – 56hz	Pain Relief level 2 – blood circulation
	P8 – 64hz	Pain Relief level 3 – blood circulation
	P9 – Sequence Loop	3 seconds P1– P8

Celsius/Fahrenheit Conversion Chart :

Celsius	Fahrenheit
10	50
20	68
30	86
40	104
50	122
60	140
70	158

FAQs :

Where is the best place on my body to use the pad?	Anywhere you feel you need TheraLOVE!
Can I use all the tech at the same time?	Yes you can!
Can I use it on my head?	This is different per technology. TENS is a NO. Low temperature InfraRed YES, PEMF YES, and red light keep to 30 minutes for anything above the neck.
Should I use this before or after my sauna session?	Used before a sauna it will open up detox pathways and help amplify the benefits of the sauna. You may even sweat more. You can also use it after. Contact: support@theraoptimalhealth.com for more info or call us at 888-416-4441.
Does it matter what heat I set it at?	Not really. Discover the one that serves you when using it before sleeping or while you sleep. Lower temps are advisable.
How long can I use the mat for?	It's different for each technology. Please see the How To Get Started page for instructions.
Who do I contact for additional support?	For product questions, email us at support@therasage.com or call 888-416-4441

Warranty

1. WHAT THIS LIMITED WARRANTY COVERS AND HOW LONG

This Product includes the product itself, and, unless specifically excluded below, all material, parts, and components related to functionality. This Limited Warranty applies only when all of these parts are genuine Therasage™ components. Cosmetic components are not covered by this Limited Warranty.

a) COVERAGE

Subject to the limitations and exclusions below, and for the duration of the applicable Warranty Period, Therasage™, and any of its affiliates, ("Seller") warrants that the Product's components will be free from material and functional non-conformities in material and workmanship under normal and reasonable consumer conditions.

b) WARRANTY PERIOD

For the one year period of this Warranty, Therasage™, at its discretion, will repair or replace any defective part, component, mechanism, or material. This warranty coverage shall include all defective parts, but shall exclude all labor charges. The Warranty Period shall commence upon purchase of the product by the end-user consumer from an authorized dealer. Each consumer is urged to register the Product through the Warranty Registration Form attached hereto. This Limited Warranty is not transferable. This Warranty can only commence once proof of purchase (or gift) is provided.

2. WHAT THIS LIMITED WARRANTY DOES NOT COVER

a) GENERALLY This Limited Warranty does not include non-conformities or damages attributable to or arising from: • Minor scratches or minor visual imperfections • Negligence; improper use, installation, finishing, maintenance operation inconsistent with Seller's recommendations and written instructions that are generally available in Seller's Product Manual as updated by bulletins or other written communications or on Seller's Worldwide website www.therasage.com • Improper storage, including inadequate shelter or inadequate venting of shipping wrap in humid locations • Misapplication of, commercial use, or use other than generally intended • Normal wear and tear • extraordinary force majeure.

b) CERTAIN CONDITIONS This Limited Warranty does not guarantee health claims. Further, this Limited Warranty does not cover damage attributable to or caused by acts of God that include, but are not limited to, power surges, stresses, fire and other conditions that exceed Product designs. Seller is not responsible for claims or damages caused by water or condensation infiltration.

3. THIS LIMITED WARRANTY'S EXCLUSIVE REMEDY

If The Product or any components fail to meet this Limited Warranty, Seller's sole obligation is to either (as Seller elects): a) repair the component (color matching not guaranteed); or b) provide replacement components to the Warranty Holder or the Seller's dealer designated (color matching not guaranteed); or c) refund the Warranty Holder's purchase price (the lesser of the original Product/ component purchase price or the original catalog list price). Repaired or replaced components are warranted only on the same terms and for the remainder of the Warranty Period. Seller reserves the right to discontinue or change any Product. If the Product or component is not available, Seller may select and provide a replacement Product or component of equal quality and price. This is the Warranty Holder 's sole and exclusive remedy for the Product under this Limited Warranty. By example but not limitation, this Limited Warranty does not cover the following costs and expenses: (i) labor; (ii)shipping/ freight expenses to return the Product to Seller, or Seller's designated service center; (iii) normal maintenance; or (iv) consequential, special, or indirect losses or damages of any kind.

4. DISCLAIMER OF WARRANTIES

THIS LIMITED WARRANTY IS IN LIEU OF AND EXCLUDES ALL OTHER WARRANTIES NOT EXPRESSLY SET FORTH HEREIN, EXPRESSED OR IMPLIED BY OPERATION OF LAW OR OTHERWISE, INCLUDING BUT NOT LIMITED TO ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. TO THE EXTENT THAT ANY IMPLIED WARRANTIES MAY NONETHELESS EXIST BY OPERATION OF LAW, SUCH WARRANTIES ARE LIMITED TO THE DURATION PROVIDED BY LAW. SOME STATES DO NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS, SO THE ABOVE LIMITATIONS MAY NOT APPLY. SELLER DOES NOT AUTHORIZE ANYONE TO CREATE FOR IT ANY OBLIGATION OR LIABILITY IN CONNECTION WITH PRODUCTS.

5. LIMITATION OF LIABILITY

SELLER'S SOLE LIABILITY UNDER THIS LIMITED WARRANTY IS REPLACEMENT, REPAIR, OR REFUND OF THE PURCHASE PRICE AS SET FORTH ABOVE. IN NO EVENT WILL SELLER BE LIABLE FOR INCIDENTAL, CONSEQUENTIAL, INDIRECT, SPECIAL, OR PUNITIVE DAMAGES INCLUDING, BUT NOT LIMITED TO, DAMAGE OF ANY KIND TO A PREMISES, LOSS OF PRODUCT USE, LABOR, LOST PROFITS/REVENUE, INTEREST, LOST GOODWILL, WORK STOPPAGE, IMPAIRMENT OF OTHER GOODS OR WORK, INCREASED OPERATING EXPENSES, EMOTIONAL DISTRESS CLAIMS, LOSS OF CONSORT CLAIMS, OR CLAIMS OF THIRD PARTIES FOR SUCH DAMAGES, WHETHER BASED ON CONTRACT, WARRANTY, TORT (INCLUDING, BUT NOT LIMITED TO, STRICT LIABILITY OR NEGLIGENCE) OR OTHERWISE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGES. SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES, SO THE ABOVE LIMITATION OR EXCLUSION MAY NOT APPLY. THIS LIMITED WARRANTY PROVIDES SPECIFIC LEGAL RIGHTS BUT WARRANTY HOLDER MAY HAVE OTHER RIGHTS WHICH VARY FROM STATE TO STATE. IF THIS LIMITED WARRANTY IS DEEMED TO HAVE FAILED OF ITS ESSENTIAL PURPOSE, IN NO EVENT WILL SELLER'S ENTIRE LIABILITY EXCEED THE LESSER OF THE PRODUCT'S OR THE NON-CONFORMING COMPONENT'S PURCHASE PRICE.

Warranty Registration

Did you purchase from us?

If you ordered your products through the Therasage Online store, your products are *automagically* registered for warranty. We all thank you for your patronage!

Did you purchase from a third party?


If purchased from a third party, in order to maximize and benefit from our Limited Warranty, don't worry we still have you covered.

Claims

Claims must be initiated during the Warranty Period. To initiate a claim, please contact Therasage™

Attn: Customer Service
10050 Remmet Ave
Chatsworth, California 91311
or by writing technicalsupport@therasage.com.

The claimant will be required to provide proof of purchase and the date of Product purchase (unless the Product was registered at the time of purchase) and may be required to return the Product or component to the Seller (at the Claimant's expense).

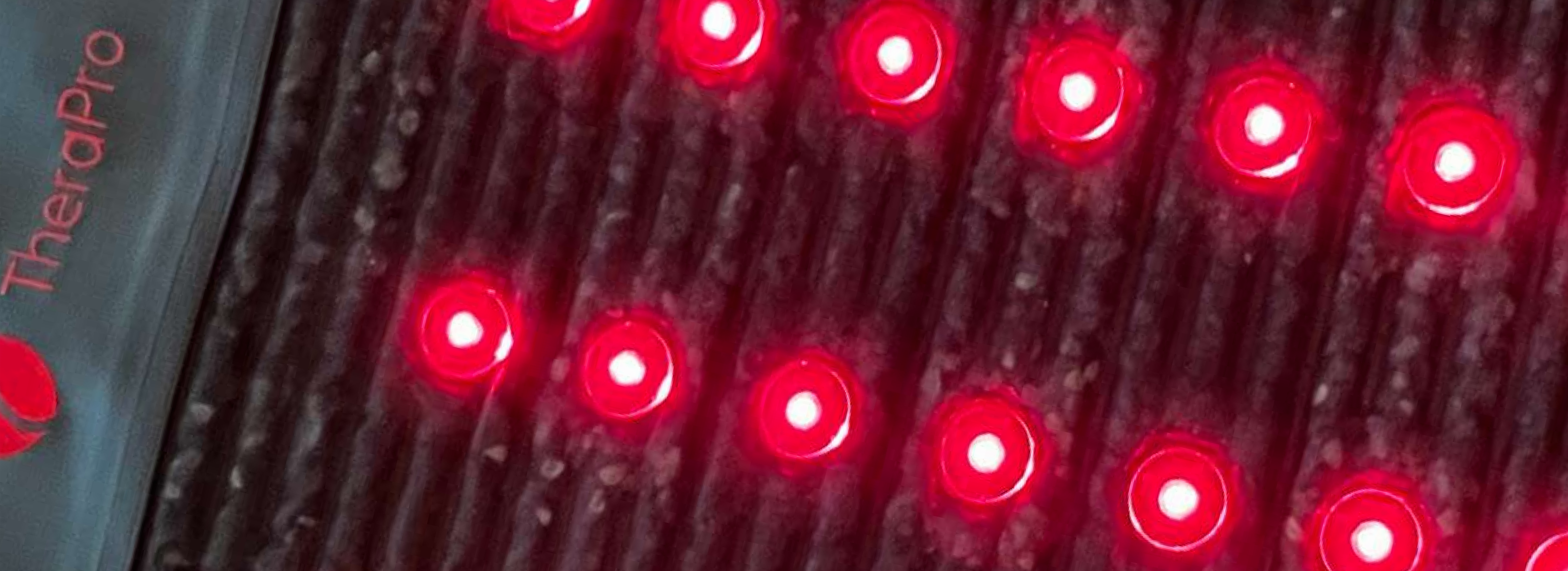

Product Related Support Email: support@therasage.com
For Health Related Questions, Email: support@theraoptimalhealth.com
Phone: 888.416.4441

Website  www.therasage.com
YouTube  Therasage
Instagram  Therasageinfrared
Facebook  Therasage
TikTok  Therasage